

# RADA EMBRYO TRANSFER CATHETER


The Embryo Transfer Catheter is a specialized and precise medical device used in the final stage of IVF/ICSI procedures for the **controlled and atraumatic transfer of embryos into the uterine cavity**.



**Precise Embryo Transfer**  
Controlled And Gentle Embryo Delivery During IVF/ICSI Procedures.



**Two-Piece Design**  
Outer Sheath And Soft Inner Catheter For Smooth Navigation.

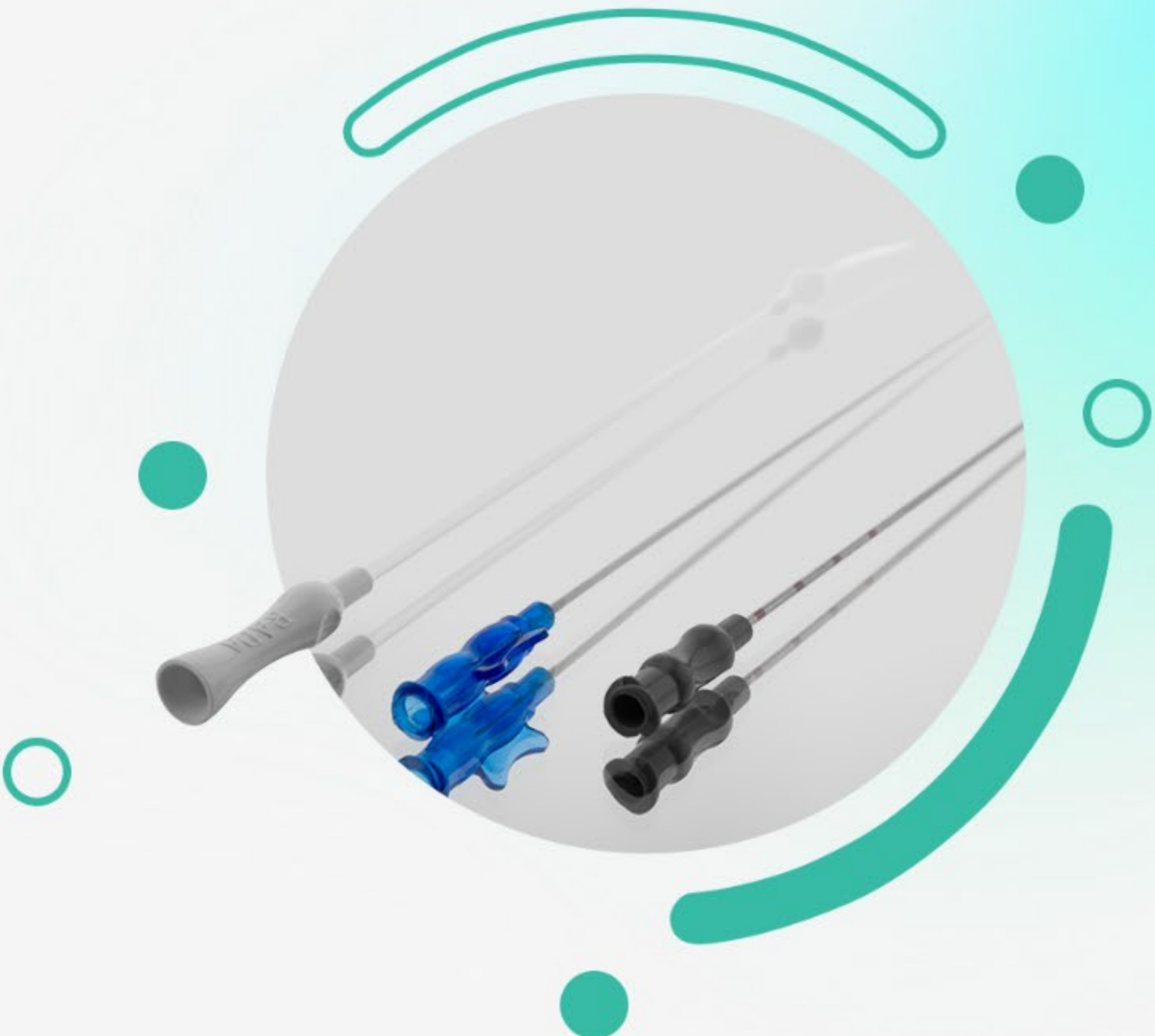


**Guided Placement**  
Accurate Catheter Positioning Under Ultrasound Guidance.



**Atraumatic Flexibility**  
Soft Tip And Adjustable Depth For Diverse Anatomies.

Featuring a **two-piece design consisting of an outer guide sheath and a soft inner transfer catheter**, this device facilitates smooth passage through the cervical canal and accurate positioning of the catheter tip under ultrasound guidance. Its atraumatic tip design, controlled flexibility, and adjustable insertion depth make it a reliable option for use across a wide range of patient anatomies and embryo transfer techniques.

## KEY FEATURES

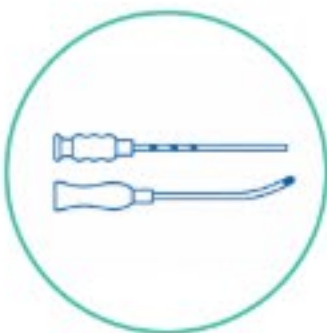

Two-piece system comprising an outer guide sheath and an inner transfer catheter



Outer guide sheath with a pre-formed curve to assist navigation through the cervical canal




Rounded, atraumatic guide tip to minimize tissue trauma and cervical irritation




Flexible outer Tube for smooth and safe cervical passage




Soft and flexible inner catheter tube for gentle and precise embryo transfer

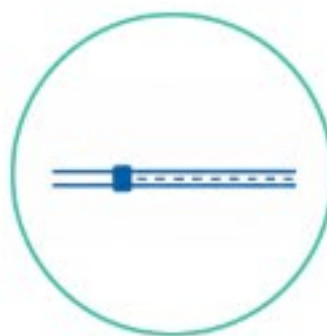

Supported inner catheter for enhanced control during insertion

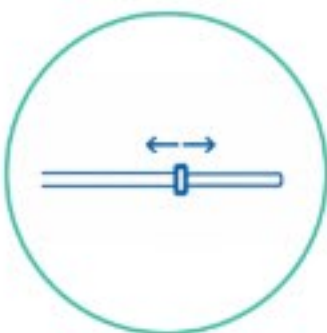

Adjustable marker for depth control and tip orientation

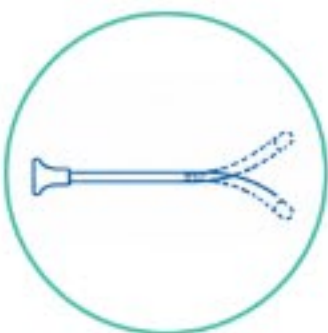

Compatible with firm mandrel, accommodating various techniques and patient anatomies



Supplied in separate sterile packaging, suitable for use in IVF and ART environments

## AVAILABLE MODELS

| IMAGE                                                                               | MODEL                  | DESCRIPTION                                                                                                         | TRANSFER TUBE LENGTH | ORDER NO. |
|-------------------------------------------------------------------------------------|------------------------|---------------------------------------------------------------------------------------------------------------------|----------------------|-----------|
|  | Simple Without Mandrel | Offers Greater Flexibility For Routine Embryo Transfer                                                              | 24 Cm                | BR344128  |
|  | Simple With Mandrel    | Provides Additional Support And Control During Insertion                                                            | 24 Cm                | BR344129  |
|  | Sono Mark              | Catheter Featuring A Metallic Ring Around The Tube For Enhanced Ultrasound Visibility                               | 24 Cm                | BR344130  |
|  | Sono Visual            | Catheter With Integrated Micro Air Bubbles Inside The Tube To Enhance Visibility Under Ultrasound Guidance          | 24 Cm                | BR344131  |
|  | Fully mark             | Catheter Featuring Metallic Rings On Both The Transfer Tube And Guide Tube For Enhanced Ultrasound Visibility       | 24 Cm                |           |
|  | Slim                   | Low-Profile Catheter With Thinner Transfer Tube And Guide Tube For Easier Insertion And Reduced Cervical Resistance | 24 Cm                |           |