

# RADA SINGLE LUMEN OOCYTE ASPIRATION NEEDLE



The Oocyte Aspiration Needle is a critical tool in the **Ovum Retrieval** process, designed for the **safe and precise collection of oocytes from follicles under ultrasound guidance**.



**Precise Oocyte Retrieval**  
Accurate Oocyte Collection Under Ultrasound Guidance.



**Enhanced Visualization**  
Echogenic Markings Improve Visibility During Guided Procedures.

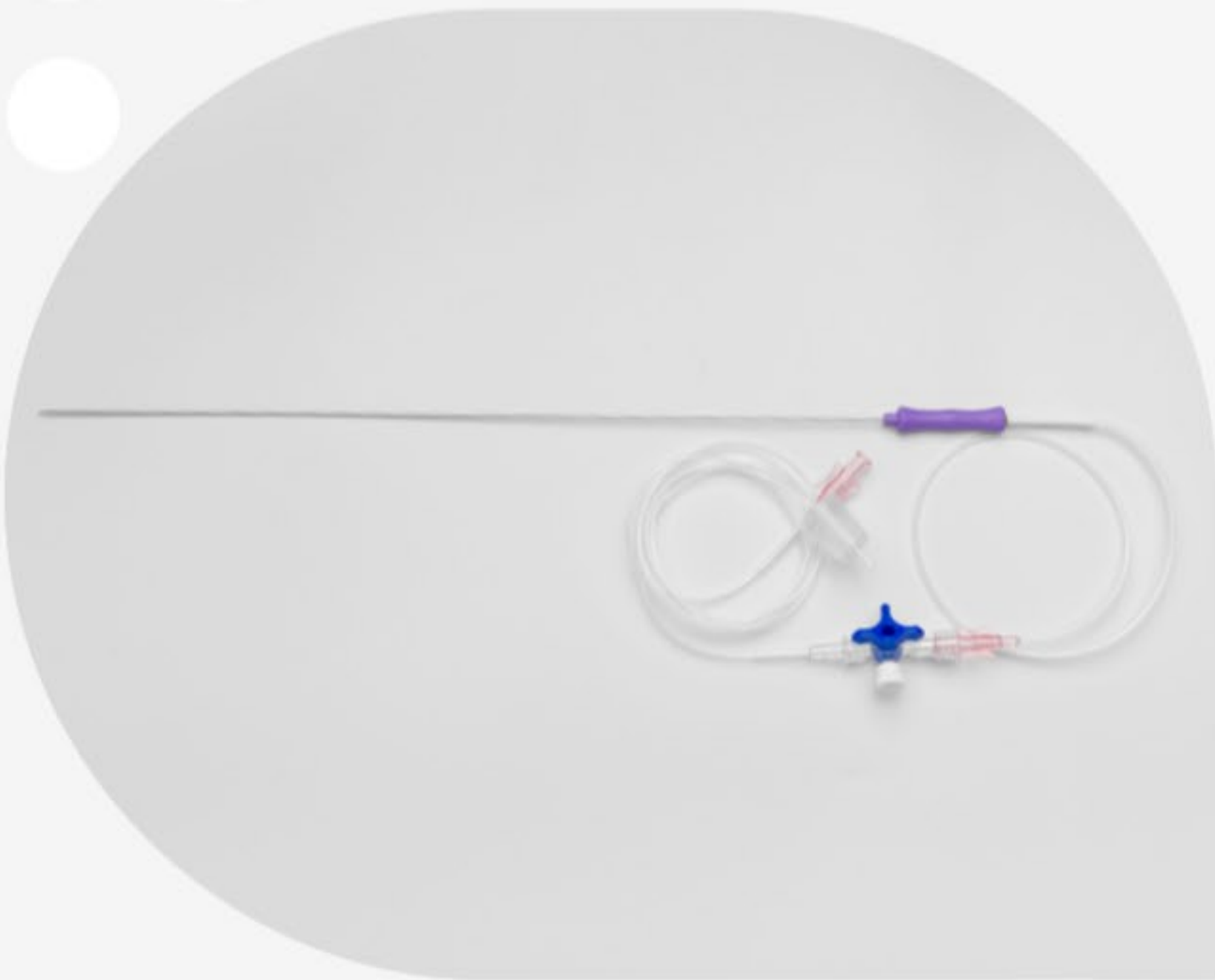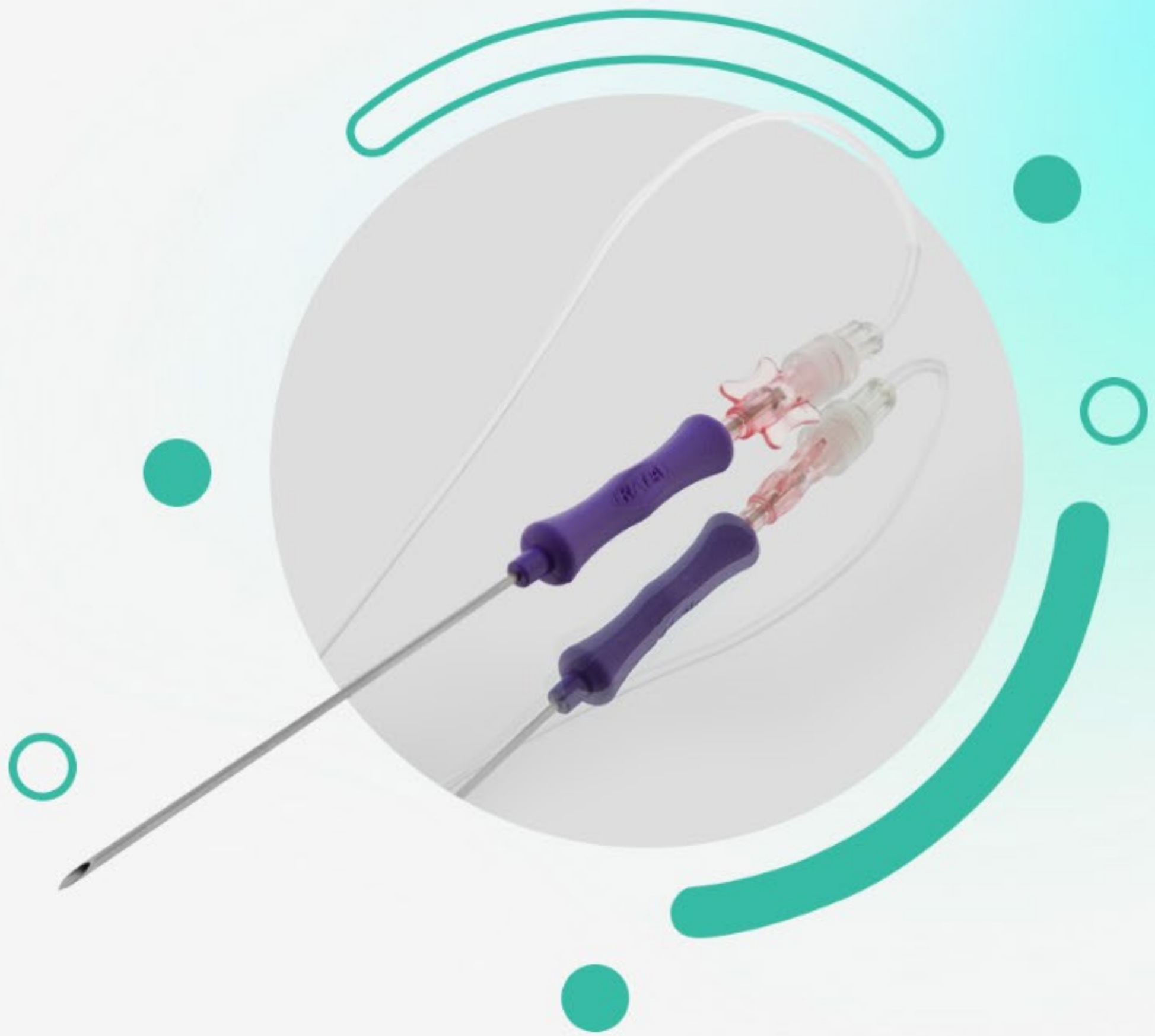


**Ergonomic Control**  
Orientation Markers Provide Confident Handling And Positioning.




**Continuous Aspiration**  
Secure Aspiration Minimizes Ovarian Trauma And Protects Oocytes.

With a sharp, atraumatic tip and **echogenic markings**, it ensures excellent visualization during ultrasound-guided procedures.  
The ergonomic handle and tip orientation markers provide full control throughout the aspiration process.  
Its single-lumen design and silicone bung connection allow for **continuous and secure aspiration**, while minimizing ovarian trauma and reducing stress on the oocytes during retrieval.


## KEY FEATURES

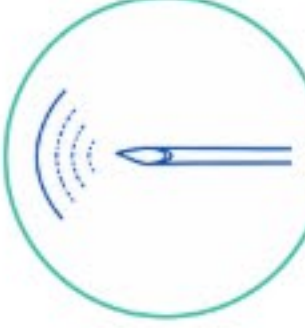


**Echogenic Tip:** Enhanced visibility under ultrasound

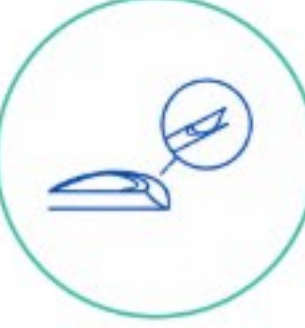


**B-bevel Tip:** Sharp, atraumatic tip for precise follicle puncture

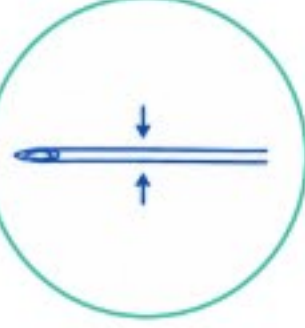


**Thin Wall Design:** Allows smaller gauge use, reducing tissue trauma and bleeding

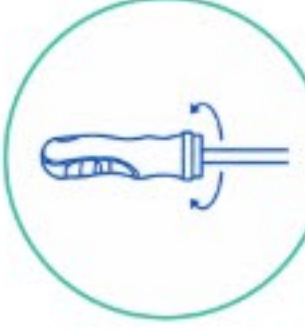

**Ergonomic Handle:** Provides control over needle rotation and tip orientation

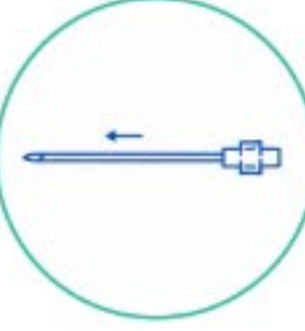

**Single Lumen:** Ensures continuous aspiration flow for minimal oocyte stress

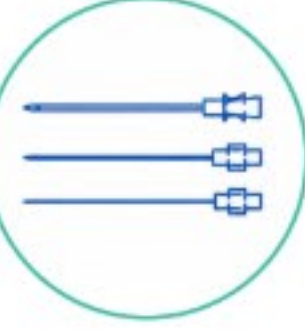

Available in multiple needle gauges and lengths to accommodate different clinical techniques and patient anatomie







**Sterile & Individually Packed:**EO-sterilized and individually packaged



**Shelf Life 3 Years:** Sterility and functionality guaranteed up to three years from manufacture

## AVAILABLE MODELS
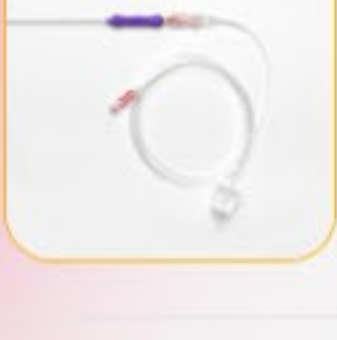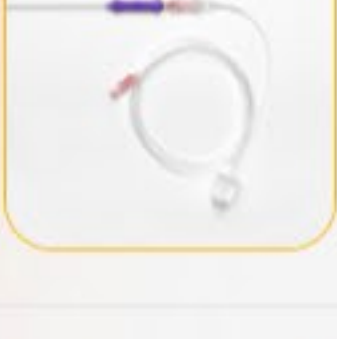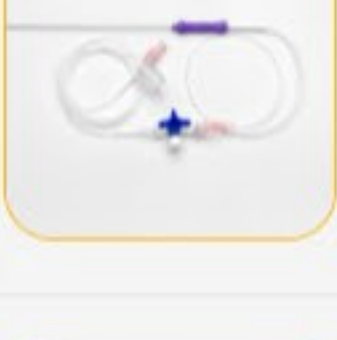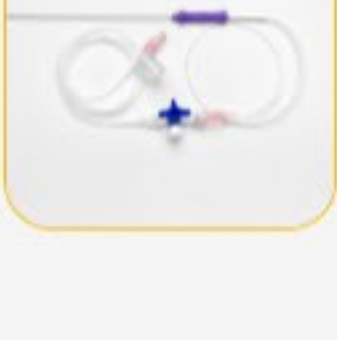
| IMAGE                                                                               | MODEL           | DESCRIPTION                                                          | ASPIRATION LENGTH | NEEDLE GAUGE | NEEDLE LENGTH | ORDER NO.                      |
|-------------------------------------------------------------------------------------|-----------------|----------------------------------------------------------------------|-------------------|--------------|---------------|--------------------------------|
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 17 G         | 35 Cm         | <a href="#">BR344314-G17</a>   |
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 17 G         | 31 Cm         | <a href="#">BR344314-S-G17</a> |
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 18 G         | 35 Cm         | <a href="#">BR344314-G18</a>   |
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 18 G         | 31 Cm         | <a href="#">BR344314-S-G18</a> |
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 19 G         | 35 Cm         | <a href="#">BR344314-G18</a>   |
|  | Simple model    | Standard needle for routine oocyte aspiration                        | 80 Cm             | 19 G         | 31 Cm         | <a href="#">BR344314-S-G19</a> |
|  | Connector model | Features a connector for easy attachment to syringes and accessories | 80 Cm             | 17 G         | 35 Cm         | <a href="#">BR344315-G17</a>   |
|  | Connector model | Features a connector for easy attachment to syringes and accessories | 80 Cm             | 17 G         | 31 Cm         | <a href="#">BR344315-S-G17</a> |
|  | 3-Ways model    | Equipped with a three-way port for easier aspiration and flushing    | 80 Cm             | 17 G         | 35 Cm         | <a href="#">BR344316-G17</a>   |
|  | 3-Ways model    | Equipped with a three-way port for easier aspiration and flushing    | 80 Cm             | 17 G         | 31 Cm         | <a href="#">BR344316-S-G17</a> |