

RADA TOWAKO METHOD




The Towako Method embryo transfer system is designed for transmyometrial embryo transfer under ultrasound guidance.



Transmyometrial Access
Direct Embryo Transfer When Transcervical Access Is Not Feasible.



Ultrasound Guided
Supports Accurate Placement Throughout The Procedure.

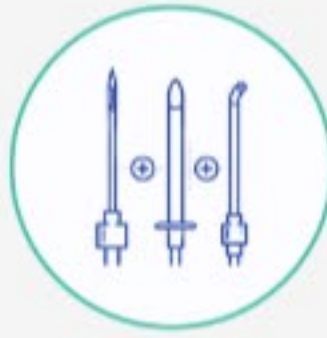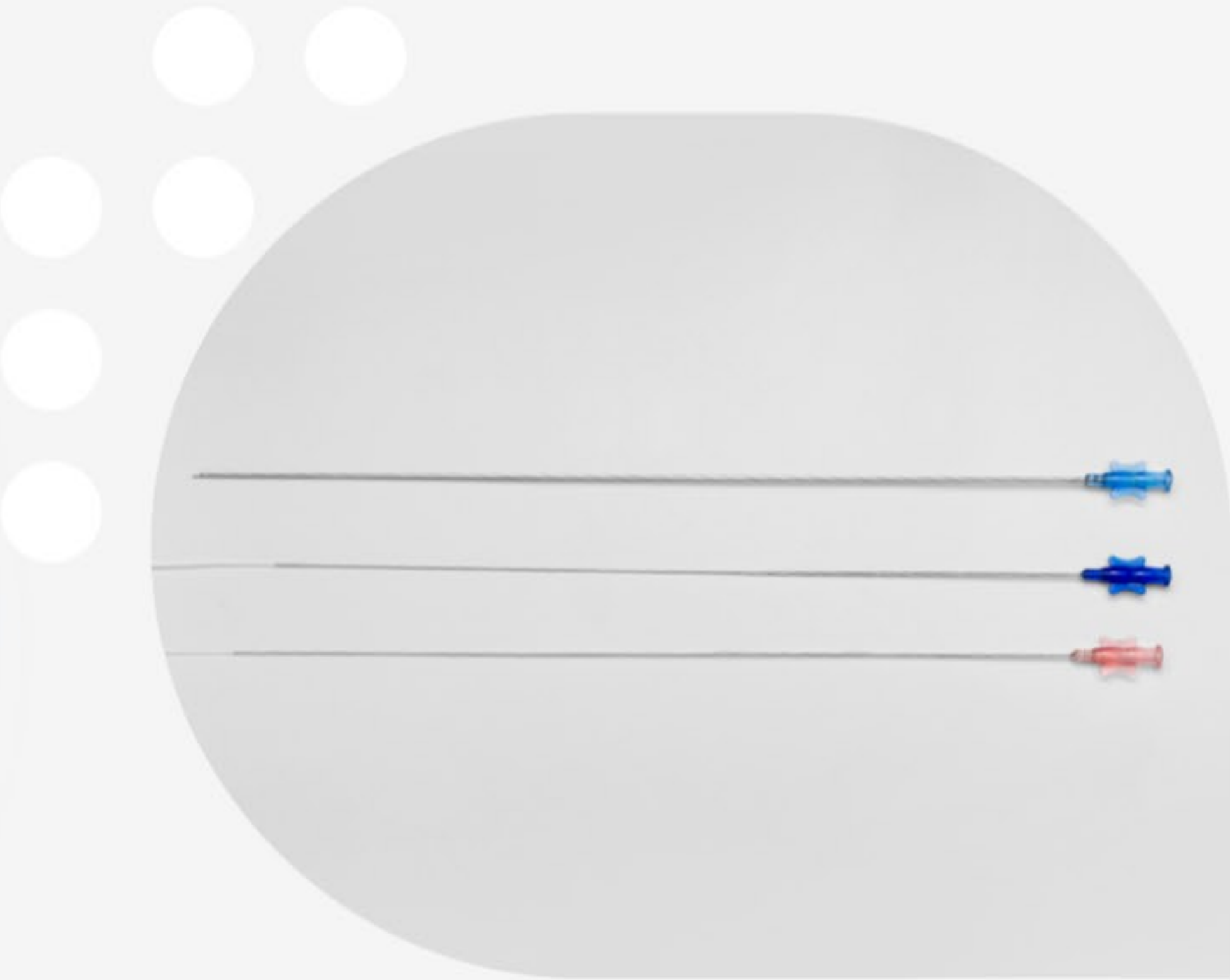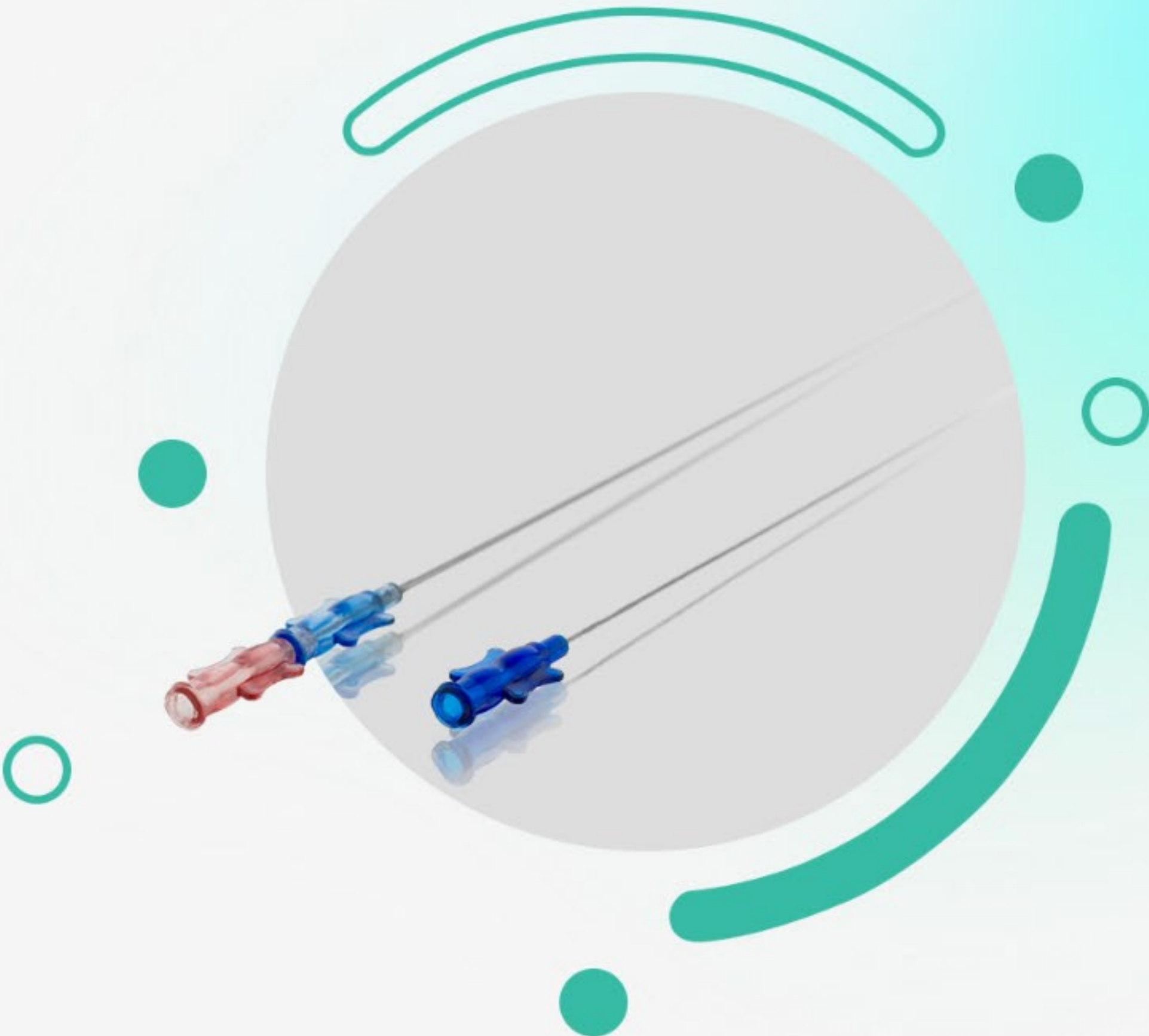

Three-Component System
Needle, Obturator, And Catheter Work Together For Controlled Access.



Precise Embryo Placement
Controlled Embryo Delivery With Reliable Positioning.

This advanced technique is indicated in cases where **conventional transcervical embryo transfer is difficult or not feasible**, providing direct access to the uterine cavity through controlled penetration of the uterine myometrium.

The system features a **three-component design**, consisting of a **Needle, obturator, and Transfer Catheter**, working together to ensure controlled access, safe penetration, and accurate embryo placement.

KEY FEATURES

Designed specifically for the Towako Method of transmyometrial embryo transfer

Intended for cases where standard transcervical embryo transfer cannot be performed

Enables direct embryo transfer through the uterine myometrium

Three-component system including:
Needle : 17G, SS304, 32 Cm
Obturator: 23G, SS304, 35 Cm
Transfer Catheter: 23G, SS304, 35 Cm

Echogenic Needle for enhanced visibility and precise control under ultrasound guidance

Smooth inner lumen catheter to support safe and atraumatic embryo transfer

Suitable for ultrasound-guided procedures

Sterile, ethylene oxide (EO) sterilized, and individually packaged

AVAILABLE MODELS

IMAGE	MODEL	DESCRIPTION	TRANSFER TUBE LENGTH	ORDER NO.
	Towako-Standard model	Three-Component System For Ultrasound-Guided Transmyometrial Embryo Transfer	Needle:17G - 32 Cm Transfer Catheter: 23G - 35 Cm	BR344166